

Med-O-Wheel Smart

for safe dosage

Dealer

Addoz Oy
 Säästöpankinkatu 2,
 30100 Forssa Finland

Tel. +358 (0)10 666 7620
 Fax +358 (0)10 296 0366

info@addoz.com
 www.addoz.com

AZ10126-0







Med-O-Wheel Smart

for safe dosage



Instructions

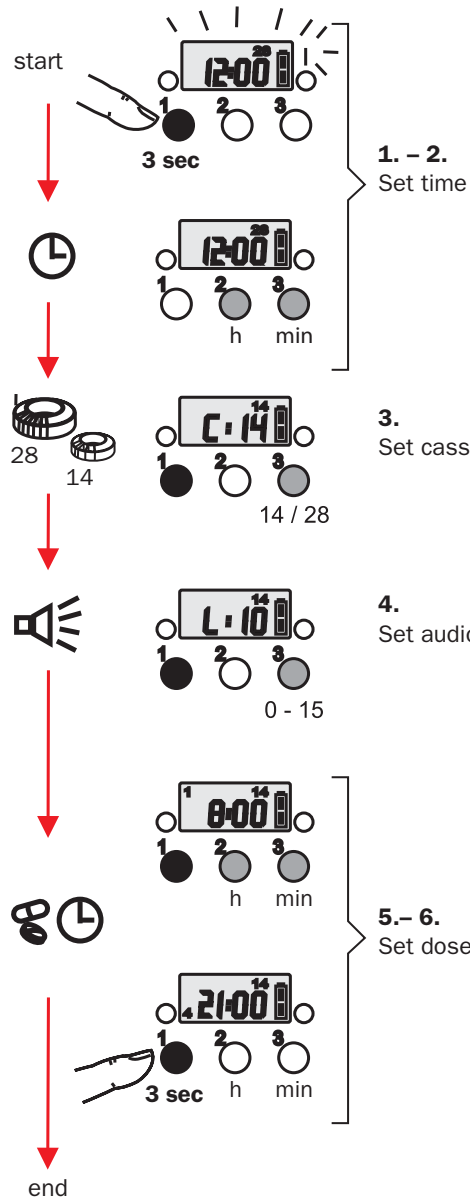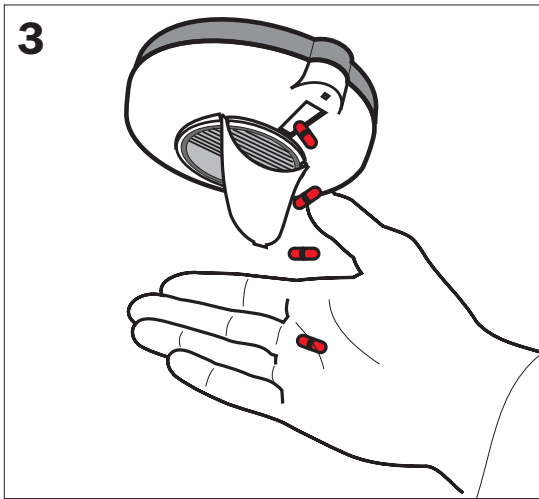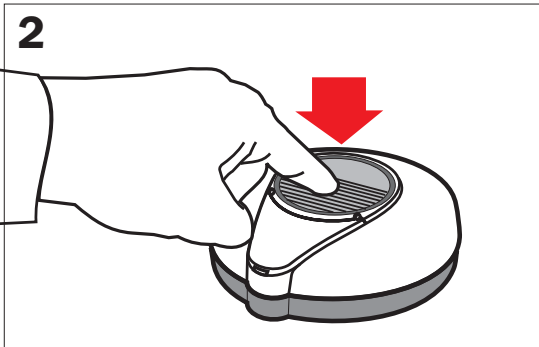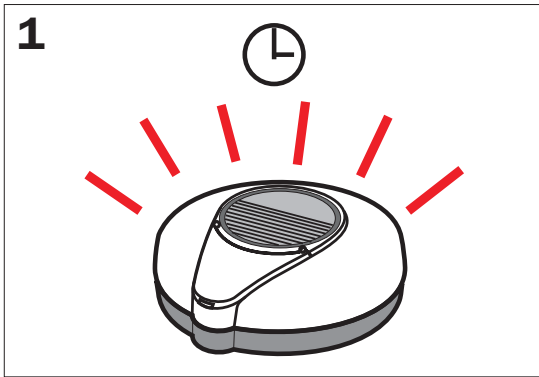

1. Press button 1 for 3 seconds - the time starts blinking.

2. Set hours with button 2 and minutes with button 3.
Accept by pressing button 1.

3. C:28/14 indicates the number of compartments.
Change by pressing button 3. Accept by pressing button 1.

4. L:15-0 indicates audio volume. Change by pressing
button 3. Accept by pressing button 1.

5. Set dose times, 1-4 times/day. Set hours with button 2
and minutes with button 3. NB! If medication is to be taken
twice a day, for example, set dose time 3 to --:-- (no time).
Accept by pressing button 1.

6. Press button 1 for 3 seconds. Cassette will automatically
return to initial position and clock is ready to use.

When the lid is closed, press on the display to check
the next dose time.